




Scan this QR code for
product information, videos,
and assembly assistance.

MEDHUB BATTERY MOUNTING BRACKET

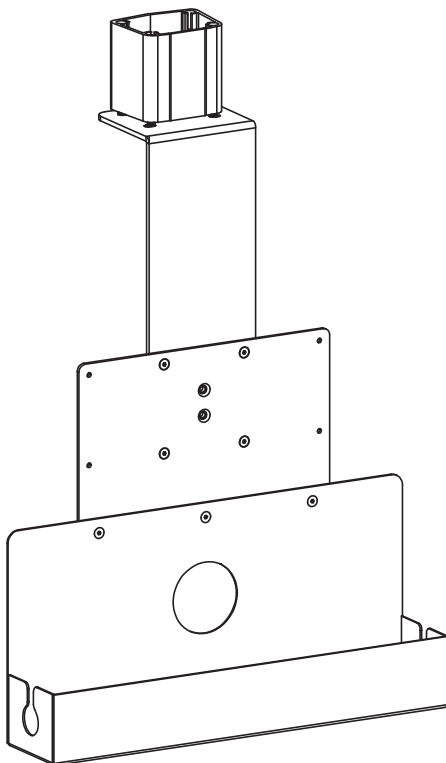

HAVE QUESTIONS?

Monday - Friday from 5:00am - 5:00pm PST



support@mount-it.com



(855) 925-5668




www.mount-it.com

For Medical Carts with Drawers

Step 1

Assemble the Mount


- Connect the L-Bracket (#C) to the Rectangular Bracket (#B) using the M4x8 Countersunk Bolts (#E) with the L-Bracket facing away from the Rectangular Bracket.
- Connect the Rectangular Bracket (#B) to the Cable Tray (#A) using the M4x8 Countersunk Bolts (#E).

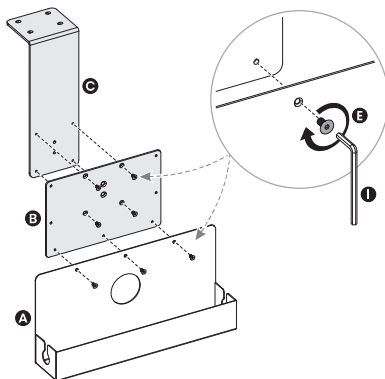

Step 2

Attach the Battery Power Unit

- Attach the power unit for the MedHub battery pack accessory (not included) to the Rectangular Bracket (#B) using the Bolts (#H).

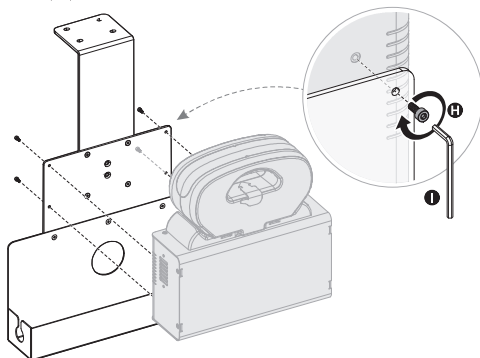


Step 3a

Mount the Bracket

For MedHub Carts with a deeper desktop tray:

- Insert the Upper Post D into the slot on the top of the tray.
- Secure the bracket from beneath the tray using the longer M6x35 Bolts (#G).

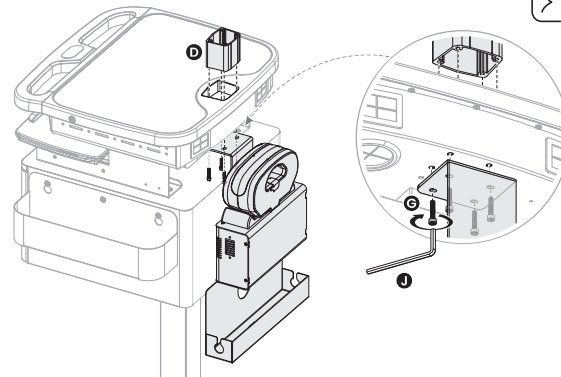


Step 3b

Mount the Bracket

For MedHub Carts with a shallow desktop tray:

- Insert the Upper Post D into the slot on the top of the tray.
- Secure the bracket from beneath the tray using the shorter M6x16 Bolts (#F).

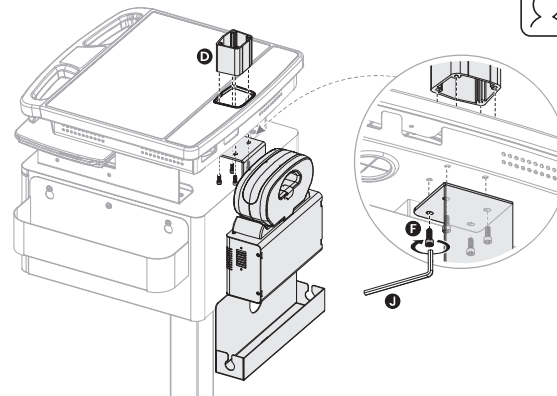


For Medical Carts without Drawers

Step 1

Assemble the Mount

- Connect the L-Bracket (#C) to the Rectangular Bracket (#B) using the M4x8 Countersunk Bolts (#E) with the L-Bracket forward from the Rectangular Bracket.
- Connect the Rectangular Bracket (#B) to the Cable Tray (#A) using the M4x8 Countersunk Bolts (#E).

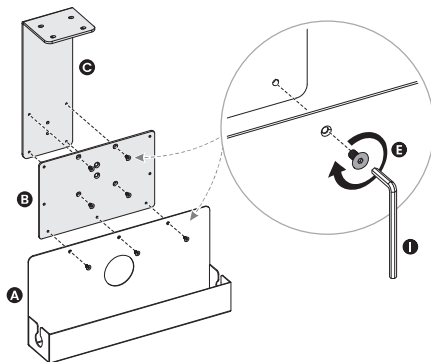


Step 2a

Mount the Bracket

For MedHub Carts with a deeper desktop tray:

- Insert the Upper Post D into the slot on the top of the tray.
- Secure the bracket from beneath the tray using the longer M6x35 Bolts (#G).

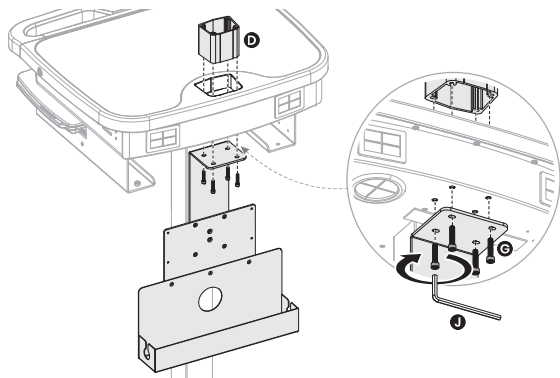


Step 2b

Mount the Bracket

For MedHub Carts with a shallow desktop tray:


- Insert the Upper Post D into the slot on the top of the tray.
- Secure the bracket from beneath the tray using the shorter M6x16 Bolts (#F).

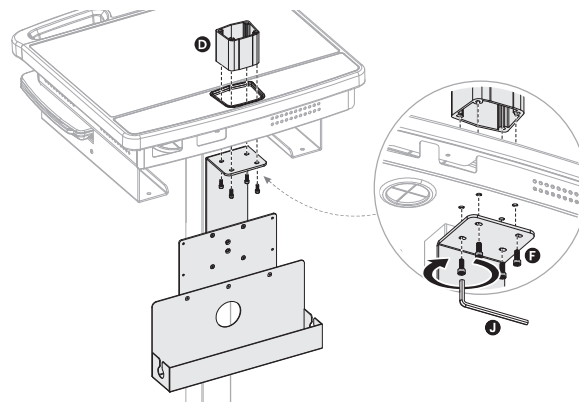

Step 3

Attach the Battery Power Unit

- Attach the power unit for the MedHub battery pack accessory (not included) to the Rectangular Bracket (#B) using the Bolts (#H).

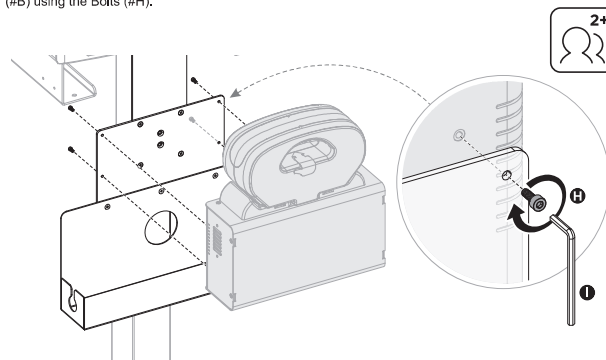

Visit us at **www.mount-it.com** for product registration and warranty information.

Our US based customer service team is standing by to answer any questions you have about your Mount-It! product.



support@mount-it.com

Available by phone **Monday - Friday from 5:00am to 5:00pm** PST.



(855) 925-5668

Chat live with an agent on our website!



www.mount-it.com