


MOUNT-IT!

HEALTH SERIES

Installation Manual

SKU: MI-16036

V2 06/26

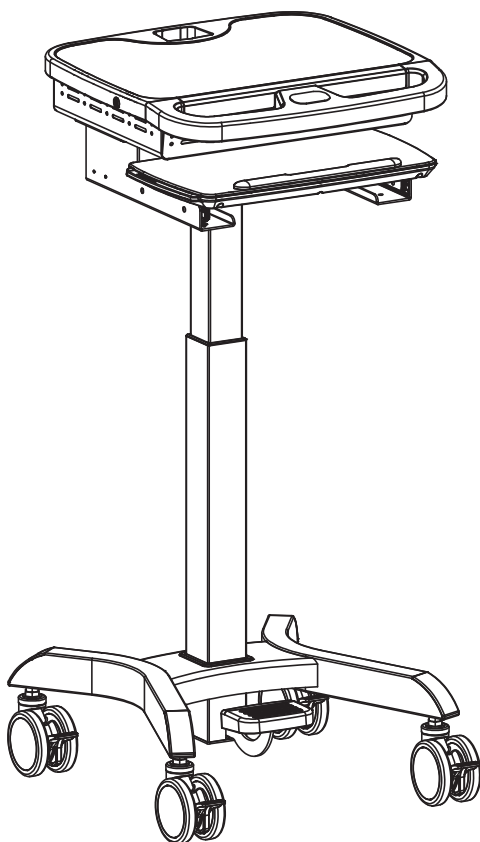

Scan this QR code for
product information, videos,
and assembly assistance.

MEDICAL CART



HAVE QUESTIONS?

Monday - Friday from 5:00am - 5:00pm PST



support@mount-it.com



(855) 925-5668




www.mount-it.com

Step 1


Attach the Casters

- Turn the Workstation (#A) upside down and attach the Casters (#B) using the Wrench (#D).

Warning: Do not hand tighten. Securely tighten both nuts on all casters and retighten after any level adjustments are made to the casters.

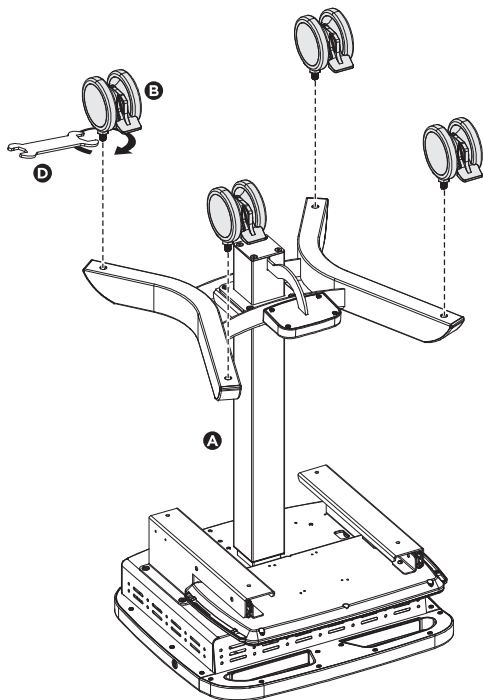


- Turn the workstation upright and raise the cart by holding each side of the platform then stepping on the foot pedal.

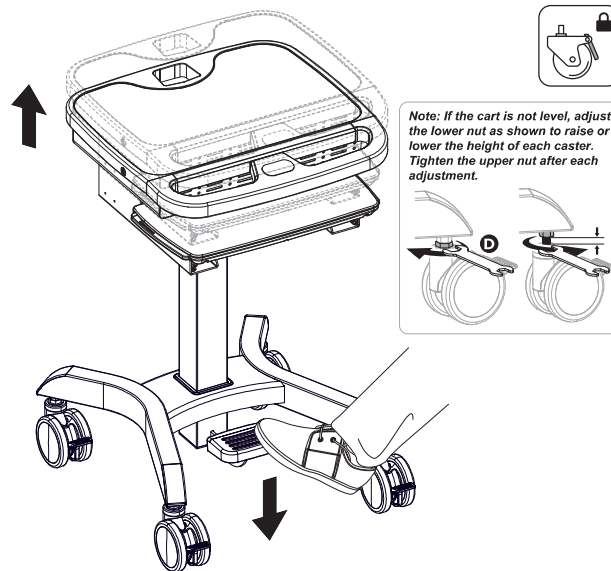

Step 2

Attach the Keyboard Stopper

- Pull the keyboard tray out. Peel the adhesive backing off the Keyboard Stopper (#C) and press it firmly into place on the end of the keyboard tray.



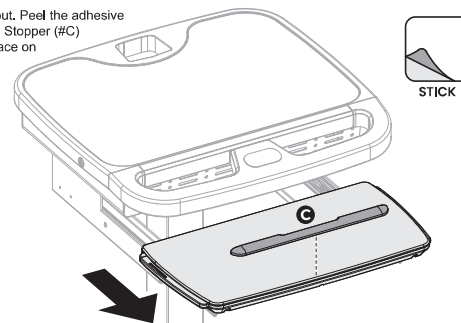


Step 3

Open the Cable Compartment

- Insert the Keys (#E) into the locks on both sides of the workstation and turn them counter-clockwise to unlock. Lift the front of the cover up, then off the workstation to access the cable compartment.


Note: The keys must be inserted into both locks at the same time. The keys will be unable to be removed from the locks while they are in the unlocked position.

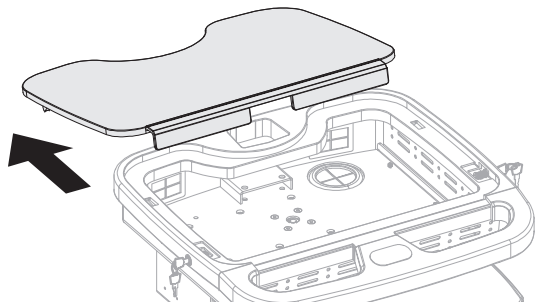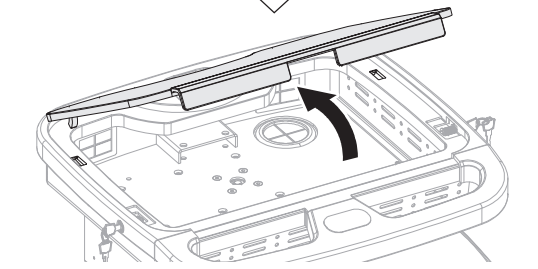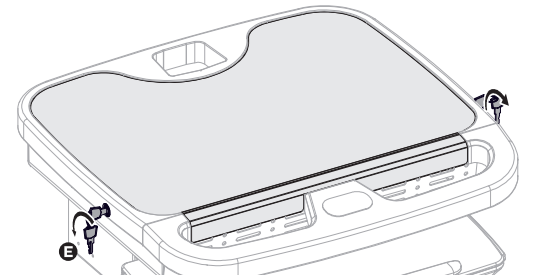

Step 4

Install the Power Strip



- Use the Adhesive Pad (#G) to secure the Power Strip (#F) into place inside the compartment.
- Route the power cable out through the port in the bottom of the compartment.

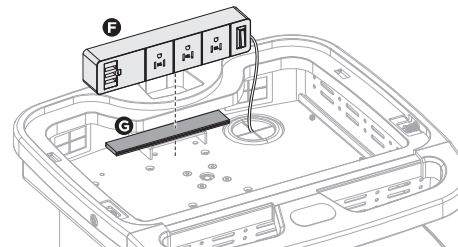


Step 5

Close the Cable Compartment

- Place the cover back onto the top of the cable compartment by first aligning the tabs with the holes near the back of the cover, then lowering it into position.
- Turn the keys clockwise to lock the compartment again.

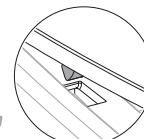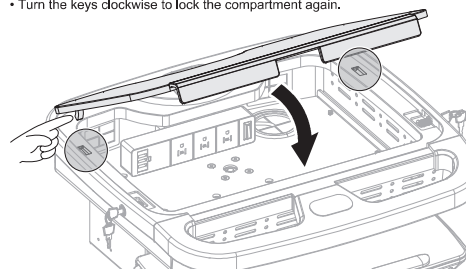

Align with the holes, and then install the desktop panel as shown below.

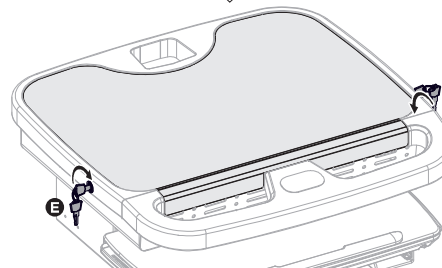

Accessories & Mounts ➡		MI-16051	MI-16052*	MI-16053	MI-16054	MI-16055	MI-16005SW	MI-16006	MI-16007	MI-16008
Carts ↓		CPU Holder	Mini CPU Holder	Scanner Holder	Wipes Holder	Printer Holder	Swivel Monitor Mount	Full Motion Monitor Arm	Full Motion Laptop Arm	Tablet Arm
MI-16036	Medhub Mobile Laptop Cart	✓	✓	✓	✓	✓	✓	✓	✓	✓
MI-16037	Medhub Mobile Laptop Cart with Drawers	□	✓	included	✓	✓	✓	✓	✓	✓
MI-16040	Medhub Motorized Mobile Laptop Cart	✓	✓	✓	✓	✓	✓	✓	✓	✓

*Requires MI-16005SW or MI-16006 monitor mount

Visit us at **www.mount-it.com** for product registration and warranty information.

Our US based customer service team is standing by to answer any questions you have about your Mount-It! product.



support@mount-it.com

Available by phone **Monday - Friday from 5:00am to 5:00pm PST.**



(855) 925-5668

Chat live with an agent on our website!



www.mount-it.com