
MOUNT-IT!
HEALTH SERIES

Installation Manual
SKU: MI-16040
V1 06/26

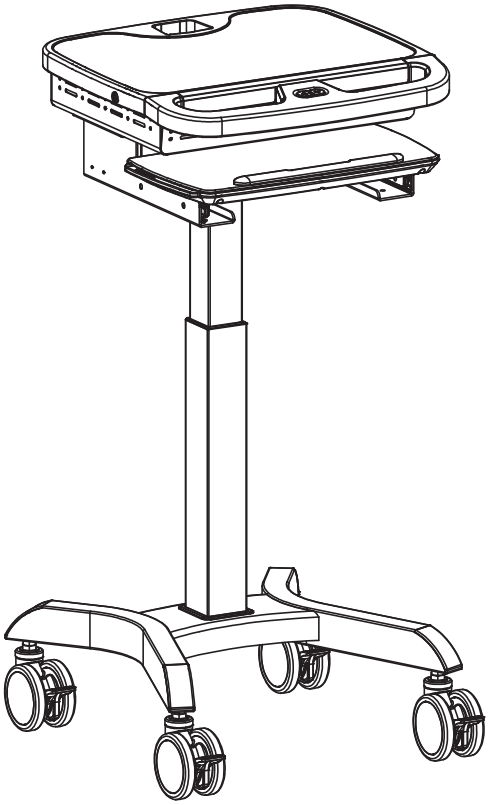Scan this QR code for product information, videos, and assembly assistance.

MEDHUB MOTORIZED COMPACT MOBILE WORKSTATION

HAVE QUESTIONS?

Monday - Friday from 5:00am - 5:00pm PST

support@mount-it.com (855) 925-5668 www.mount-it.com

Thank you for choosing Mount-It!

EN **IMPORTANT:**
Failure to read, thoroughly understand, and follow all instructions can result in serious personal injury, damage to equipment, or voiding of factory warranty.

INSTRUCTIONS OF SAFETY AND WARNINGS:
- Make sure that the mounting surface is strong enough to handle a mounted product and equipment.
- DO NOT EXCEED THE MAXIMUM LISTED WEIGHT CAPACITY.
- Always use an assistant or a mechanical lifting device to safely lift and position heavy products.
- Make sure screws firmly, but do not over tighten. Over tightening can cause product damage that greatly reduces the holding power.
- Keep a clear area and distance from moving parts when using the product.
- Use only as intended. Never stand, hang on or climb on the product.
- This product is intended for indoor use only.
- Any attempts to reconstruct the construction aren't allowed.
- This product may contain small items that could be a choking hazard if swallowed. Keep children away unless the product is intended for their use and all instructions and manuals are fully reviewed and understood by them.
- Check that the product is secure and safe to use at regular intervals (at least every three months).

Should you ever have any questions or feel you require help please contact Mount-It!

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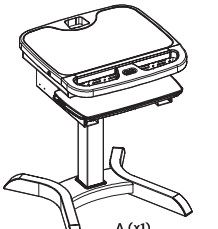
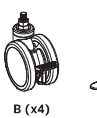
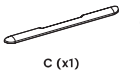
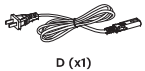
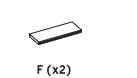
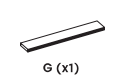
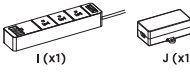
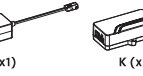
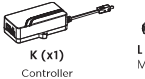
K (x1)
Controller


L (x2)
M4x8

Check carefully to make sure there are no missing or defective parts. Some hardware provided may not be used. Use the correct mounting hardware suitable for your installation.

CHOKING HAZARD. Package includes small parts.
Not for children under 3 years. Adult supervision is required.

Missing Parts? Call us at (855) 925-5668 and we will ship replacement parts free of charge.


DIMENSIONS & WEIGHT CAPACITY

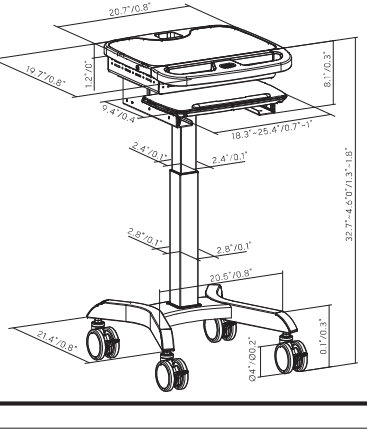Rated weight capacity is 110lbs/50kgs for the medical cart.
Caution:
Ensure all mounting surfaces are capable of holding the weight of your mount.
DO NOT EXCEED MAXIMUM LISTED WEIGHT CAPACITY. SERIOUS INJURY OR PROPERTY DAMAGE MAY OCCUR!


Step 1
Install the Casters

• Turn the Workstation (#A) upside down to access the base.
• Insert each Caster (#B) into the base and tighten securely using the 19–21mm Wrench (#H).
Warning: Do not hand tighten. Securely tighten both nuts on all casters and retighten after any level adjustments are made to the casters.

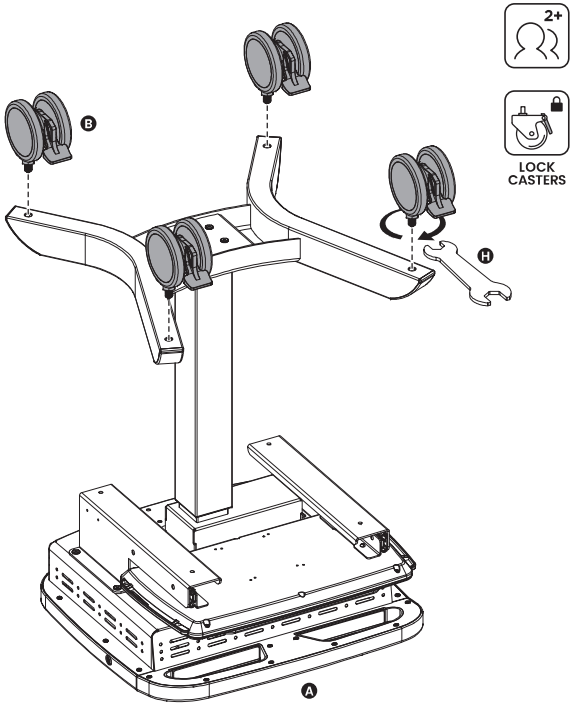
Step 2
Remove the Top Tray


• Turn the Workstation (#A) upright.
• Use the Keys (#E) to unlock both sides of the tray, then lift the front of the tray and remove it from the Workstation (#A).

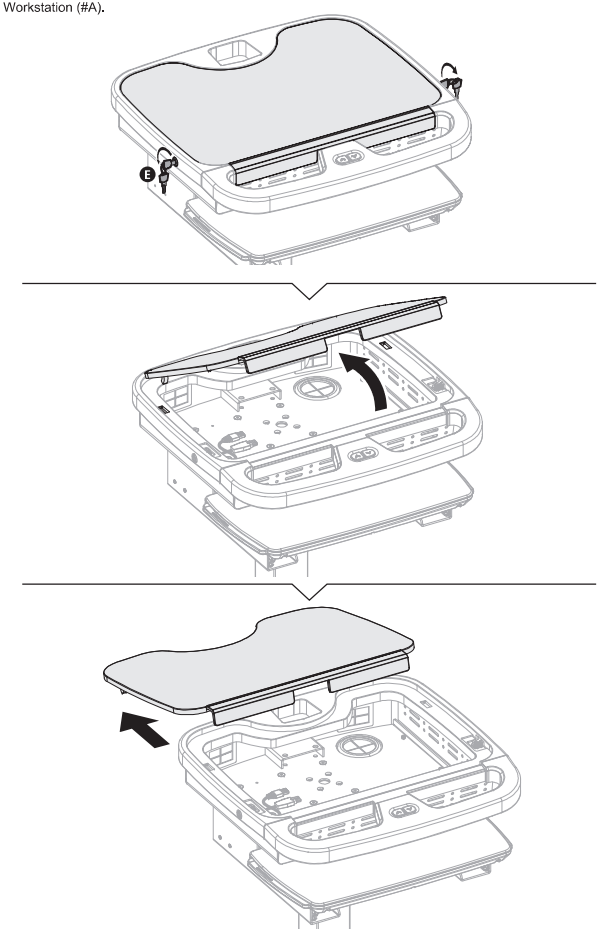
Step 3
Install Power Components

• Install the AC Adapter (#J) using the Bolts (#L) and the Small Adhesive Pad (#F).
• Install the Controller (#K) using the remaining Small Adhesive Pad (#F).

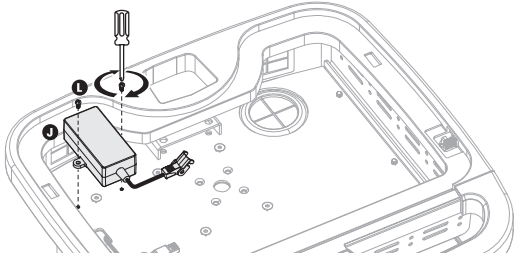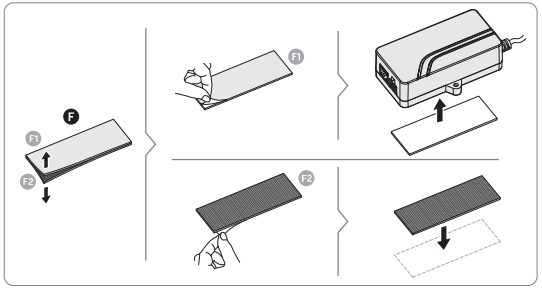
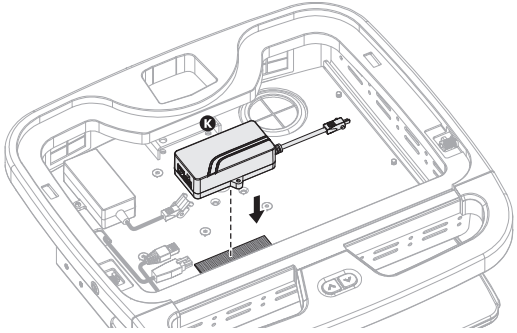
Step 4
Connect Internal Cables


• Inside the workstation tray, connect the cables as shown in the diagram, ensuring all connections are fully seated.

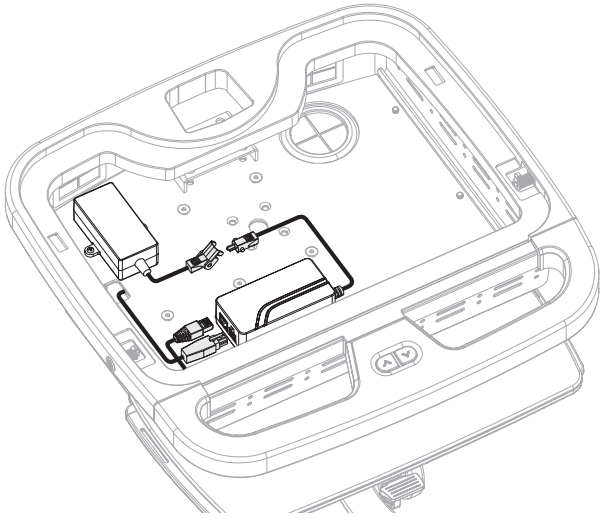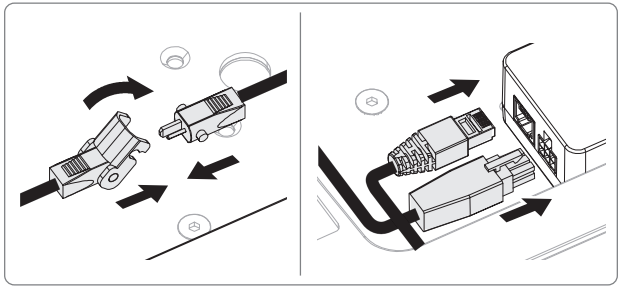
Step 5
Install Power Components


• Insert the Power Cable (#D) through the rear slot of the Workstation (#A) and connect it to the internal power adapter.
• Peel and apply the Large Adhesive Pad (#G) as shown, secure the Power Strip (#I) on top, then route its cable through the bottom opening of the tray.

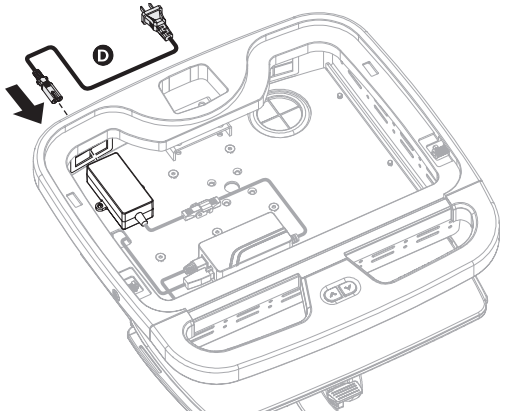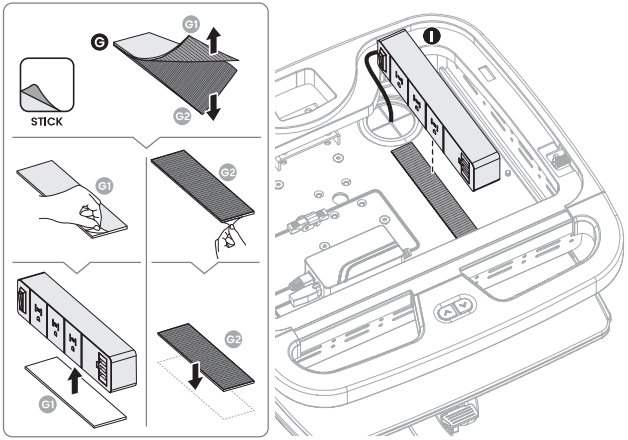
Step 6
Reinstall the Top Tray

• Place the tray back onto the Workstation (#A), aligning it with the mounting holes as shown, then use the Keys (#E) to lock both sides.

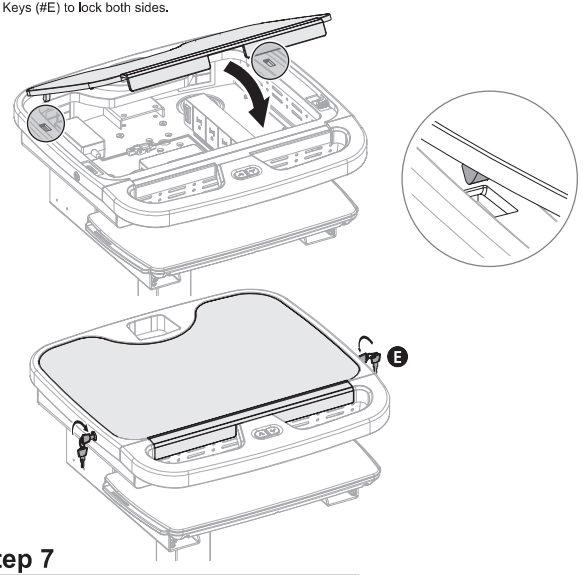Step 7
Install the Keyboard Stopper


• Pull out the keyboard tray, then peel and apply the Keyboard Stopper (#C) to the front edge as shown.

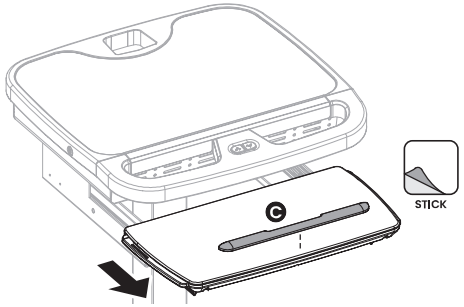
Step 8

Leveling and Use

- The mouse tray can be extended from either side of the keyboard tray.
- Lock the Casters when not in use or when a more stable workstation is needed.

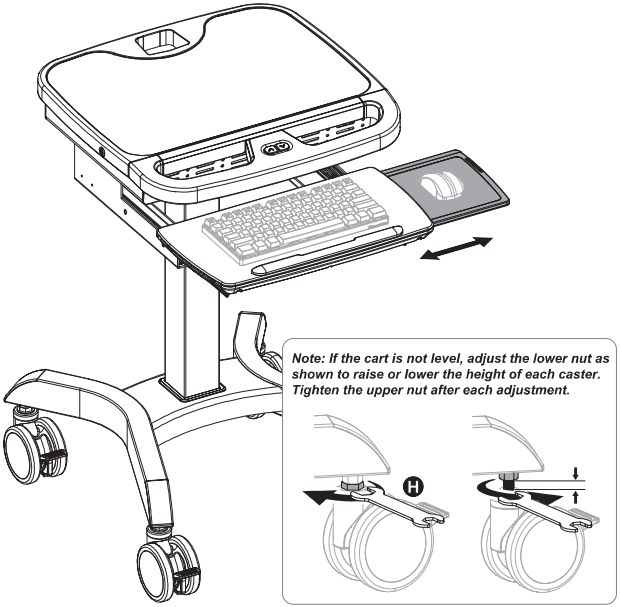Note: If the cart is not level, adjust the lower nut as shown to raise or lower the height of each caster. Tighten the upper nut after each adjustment.

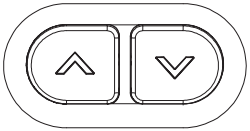1. Press the " ^ " button to raise the desk, press the " v " button to lower the desk.
2. At the lowest position, simultaneously hold the " ^ " and " v " buttons to perform the reset.

The system may be overheating when the unit beeps and can not go up or down. Stop operating and allow your desk to remain idle for approximately 18 minutes.

Accessories & Mounts →		MI-16051	MI-16052*	MI-16053	MI-16054	MI-16055	MI-16005SW	MI-16006	MI-16007	MI-16008
Carts ↓		CPU Holder	Mini CPU Holder	Scanner Holder	Wipes Holder	Printer Holder	Swivel Monitor Mount	Full Motion Monitor Arm	Full Motion Laptop Arm	Tablet Arm
MI-16036	Medhub Mobile Laptop Cart	✓	✓	✓	✓	✓	✓	✓	✓	✓
MI-16037	Medhub Mobile Laptop Cart with Drawers	☐	✓	Included	✓	✓	✓	✓	✓	✓
MI-16040	Medhub Motorized Mobile Laptop Cart	✓	✓	✓	✓	✓	✓	✓	✓	✓

*Requires MI-16005SW or MI-16006 monitor mount

IMPORTANT:

The product is powered by electricity. In order to avoid burns, fire and electric shock, please read the instructions carefully. Failure to read, thoroughly understand, and follow all instructions can result in serious personal injury, damage to equipment, or voiding of factory warranty.

Other Cautions

- Don't clean the product with water while the power is still connected.
- Don't disassemble or replace components while the power is still connected.
- Never operate the system with a damaged cord or plug.
- Never operate the system if it's in moist environment or its electrical components have contact with liquids.
- If the control box makes unusual noise or smells, switch off the mains voltage immediately.

Disposal: The product with this mark indicates that this product shouldn't be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the local recycle service or contact your point of purchase to withdraw the product.

Visit us at **www.mount-it.com** for product registration and warranty information.

Our US based customer service team is standing by to answer any questions you have about your Mount-It! product.

Available by phone **Monday - Friday from 5:00am to 5:00pm** PST.

Chat live with an agent on our website!

support@mount-it.com

(855) 925-5668

www.mount-it.com